

5280  Creations  
LAB SPECIALISTS

# Implant Restorative Options

## Cement retained Custom abutment:


- Crown and custom abutment design NO screw access hole in crown.
- This is ideal for tight interproximal clearance on smaller implant platforms, Less than perfect screw access entry on the facial or incisal edge of anterior teeth.



---

## Custom Screw-mentable Design:

- Crown with screw access hole and custom abutment will come to you in 2 pieces, for in office cement.
- This is a good option for cases that do not allow for parallel contacts/draw.

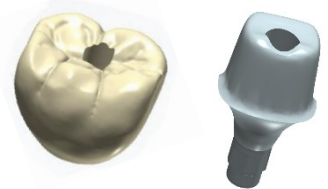

---

## Custom Screw retained Design:

- Crown and custom abutment cemented together in lab and delivered as one piece.
- This is ideal for cases that have a parallel draw or willing to adjust adjacent contacts to seat properly.



---

## Custom Cast abutment (UCLA):

- True all in one design. Porcelain is stacked directly to the anatomically casted Gold/High Nobel abutment.

OR

- A similar option we have for this method would be a CAD designed high-quality cobalt-chromium alloy for venerable screw retained restorations. (Straumann Specific)



---

## 5280ALLin1:

- Crown with screw access hole on a stock Ti base, cemented together.
- This is ideal for pmma temps for tissue training or a more economical solution for your final implant restoration.
- Posteriors only



---

**\*For further Implant case planning and extensive material options**

**contact Megan 303-280-5280 / [Megan@5280creations.com](mailto:Megan@5280creations.com)\***